# TEMPEST: Porewater SO4/Cl

# Apr-Dec (No May or June) 2024 Samples

## 2025-10-20

# Contents

0.1	Run Information	2
0.2	Pull in active porewater tracking inventory sheet from Google Drive:	2
0.3	Assess Standard Curves	3
0.4	Assess Check Standards	4
0.5	Assess Blanks	5
0.6	Assess Duplicates	7
0.7	Calculate mmol/L concentrations $\&$ salinity, add dilutions	8
0.8	Assess Analytical Spikes	9
0.9	Check if samples within the range of the standard curve $\ldots \ldots \ldots \ldots \ldots$	9
0.10	Check to see if samples run match metadata & merge info	10
0.11	Visualize Data by Plot	11
0.12	Export Processed Data	11

#### 0.1 Run Information

```
###### Run information - PLEASE CHANGE
  Sample_Year = "2024"
  Date_Run = "2025-08-13" #Date that instrument was run
  Run_by = "Zoe Read" #Instrument user
  Script_run_by = "Zoe Read" #Code user
  run_notes = "Check standard 1 is low for S04, so <80% of Sulfate Check
  Standards are within range of expected concentration. The peaks look uniform,
  so I'm not sure why the standard is low.
  " #any notes from the run
##Sample data that was entered incorrectly
# The Old ID is the original, incorrecly-entered ID and the New ID is the correct ID to change it to.
Old_ID_1 = "TMP_SW_F6_20241111_15cm"
New_ID_1 = "TMP_SW_F6_20241104_15cm"
###### File Names - PLEASE CHANGE
#file path and name for raw summary data file
 raw_file_name_cl = "Raw Data/COMPASS_TEMPEST_202404-202412_Cl.txt"
 raw_file_name_so4 = "Raw Data/COMPASS_TEMPEST_202404-202412_S04.txt"
#file path and name of processed data file
 processed_file_name = "Processed Data/COMPASS_TEMPEST_Processed_C1_S04_202404-202412.csv"
###### Log Files - PLEASE CHECK
#qaqc log file path for this year copied over from COMPASS GitHub
 Log_path = "Raw Data/COMPASS_Synoptic_C1_S04_QAQClog_2024.csv"
```


##Set Up Code - constants and QAQC cutoffs

#### 0.2 Pull in active porewater tracking inventory sheet from Google Drive:


```
##Create similar sample IDs to match with run samples
##Import Sample Data
##Fix Sample IDs entered wrong
```

### 0.3 Assess Standard Curves

## Chloride Std Curve

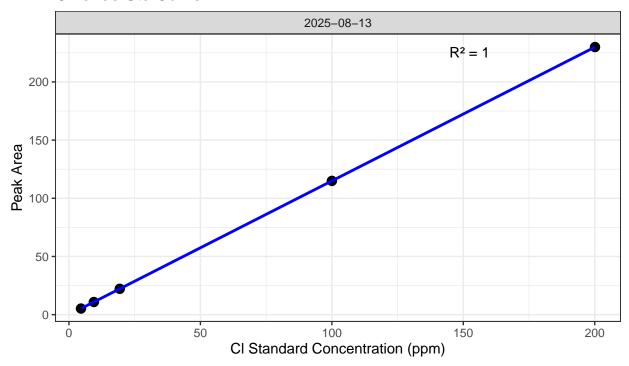


# Sulfate Std Curve

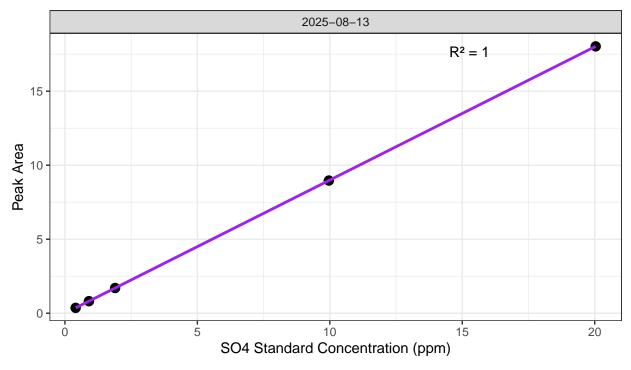

 $\mbox{\tt \#\#}$  [1] "QAQC log file exists and has been read into the code."

### Slope Drift Assessment

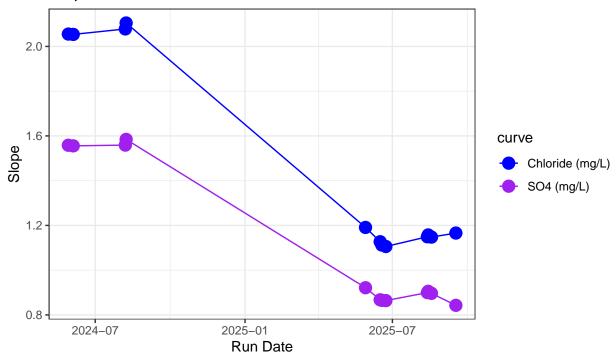


## [1] "Cl Curve r2 GOOD"

## [1] "SO4 Curve r2 GOOD"

#### 0.4 Assess Check Standards

```
## # A tibble: 5 x 5
     sample ID mean Cl sd Cl cv Cl flag Cl
##
##
     <chr>>
                 <dbl> <dbl> <dbl> <chr>
## 1 Standard 1
                  4.67 0.126 0.0270 Chloride Check Standard RSD within Range - P~
## 2 Standard 2
                  9.60 0.153 0.0159 Chloride Check Standard RSD within Range - P~
                19.5 0.296 0.0152 Chloride Check Standard RSD within Range - P~
## 3 Standard 3
                       0.188 0.00188 Chloride Check Standard RSD within Range - P~
## 4 Standard 4 100.
## 5 Standard 5 201.
                       0.852 0.00424 Chloride Check Standard RSD within Range - P~
## # A tibble: 5 x 5
     sample_ID mean_S04 sd_S04 cv_S04 flag_S04
##
     <chr>>
                  <dbl>
                          <dbl>
                                  <dbl> <chr>
## 1 Standard 1
                  0.409 0.00380 0.00929 Sulfate Check Standard RSD within Range -~
                  0.910 0.00386 0.00424 Sulfate Check Standard RSD within Range -~
## 2 Standard 2
                  1.90 0.00779 0.00409 Sulfate Check Standard RSD within Range -~
## 3 Standard 3
                  9.96 0.0194 0.00195 Sulfate Check Standard RSD within Range -~
## 4 Standard 4
## 5 Standard 5
                 20.1
                        0.0776 0.00386 Sulfate Check Standard RSD within Range -~
## [1] ">80% of Chloride Check Standards have RSD within range - PROCEED"
```


## [1] ">80% of Sulfate Check Standards have RSD within range - PROCEED"

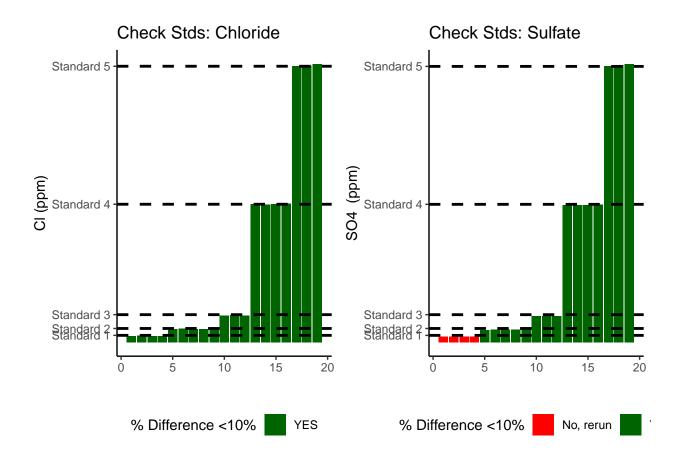

## [1] ">80% of Chloride Check Standards are within range of expected concentration - PROCEED"
## [1] "<80% of Sulfate Check Standards are within range of expected concentration - REASSESS"</pre>


#### 0.5 Assess Blanks

- ## [1] ">80% of Chloride Blank concentrations are lower 25% quartile of samples"
- ## [1] ">80% of Sulfate Blank concentrations are lower 25% quartile of samples"

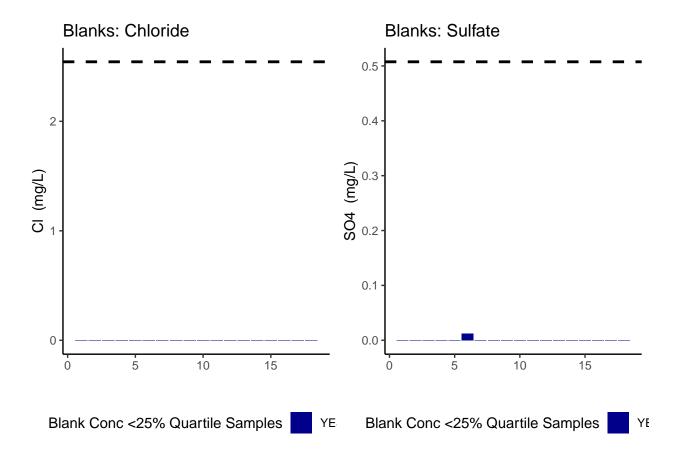


- ## Chloride blanks mean ppm:
- **##** [1] 0
- ## Sulfate blanks mean ppm:
- ## [1] 0.0006722222

## 0.6 Assess Duplicates

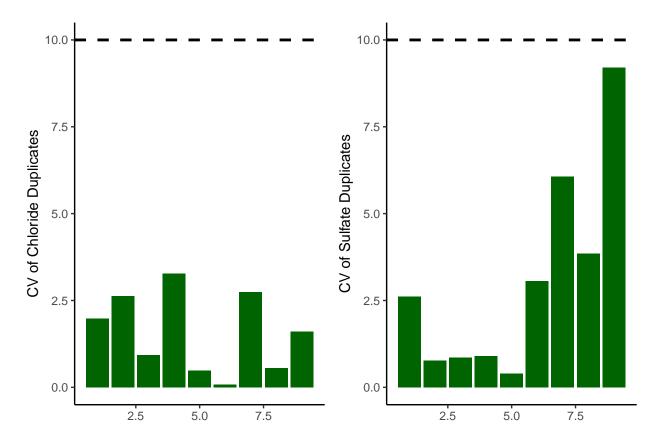

- ## [1] ">80% of Chloride Duplicates have a CV <10% PROCEED"
- ## [1] ">80% of Sulfate Duplicates have a CV <10% PROCEED"

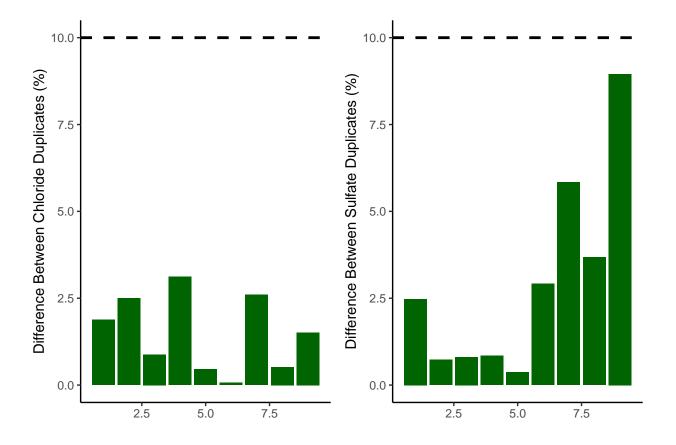


- ## [1] ">80% of Chloride Duplicates have a percent difference <10% PROCEED"
  ## [1] ">80% of Sulfate Duplicates have a percent difference <10% PROCEED"</pre>
- 0.7 Calculate mmol/L concentrations & salinity, add dilutions

```
# Convert ppm to mmol/L
all_dat$S04_Conc_mM <- (all_dat$S04_ppm / s_mw)
all_dat$Cl_Conc_mM <- (all_dat$Cl_ppm / cl_mw)

# Calculate Salinity
# calculated using the Knudsen equation
# Salinity = 0.03 + 1.8050 * Chlorinity
# Ref: A Practical Handbook of Seawater Analysis by Strickland & Parsons (P. 11)
# =((1.807*Cl_ppm)+0.026)/1000
all_dat$salinity <- ((1.8070 * all_dat$Cl_ppm) + 0.026) / 1000

#Need to determine dilution factors for your samples
#for TEMPEST this depends on the sample so...
all_dat$Dilution <- 1
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "TMP"), 50, all_dat$Dilution)
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "EST_SourceWater"), 100, all_dat$Dilution)
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "SW_SourceWater"), 100, all_dat$Dilution)</pre>
```

## 0.8 Assess Analytical Spikes

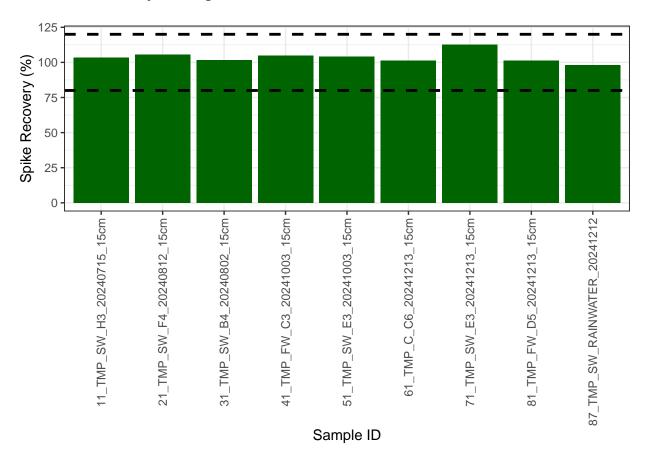

## [1] ">80% of SO4 spikes have a recovery between the high and low cutoff - PROCEED"

## 0.9 Check if samples within the range of the standard curve

## Sample Flagging

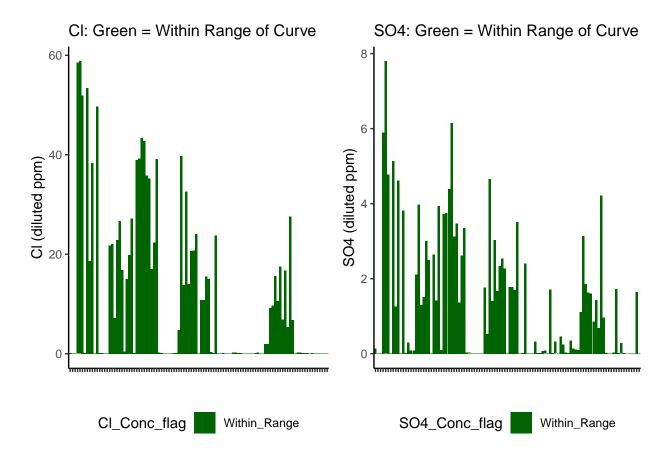


Table 1: SO4 samples

SO4_Conc_flag	Percent_samples
Within_Range	100


Table 2: Cl samples

Cl_Conc_flag	Percent_samples
Within_Range	100

## 0.10 $\,$ Check to see if samples run match metadata & merge info

## All sample IDs are present in metadata.

## 0.11 Visualize Data by Plot

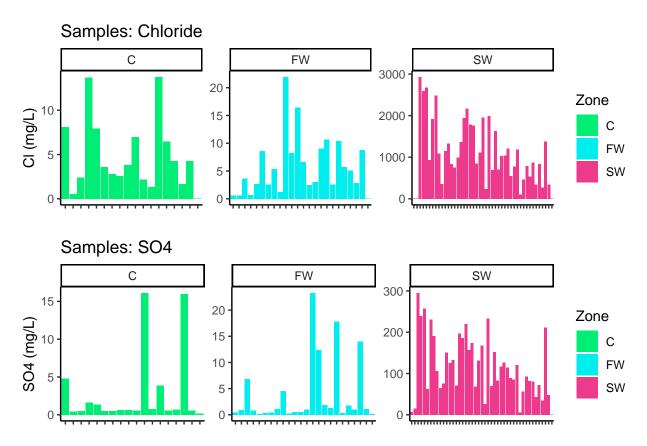

## 0.12 Export Processed Data

#end